BACtrack Go

Keychain Breathalyzer BT-KC20 — User Guide

Select a language

EN English page 2
DE Deutsch . . . page 12
CS Čeština . . . page 22
SK Slovenčina page 32
PL Polski page 42
HU Magyar . . . page 52
RO Română . . . page 62
SL Slovenščina page 72
BG Bulgarian . page 82

BACtrack® Go

Keychain Breathalyzer

Model BT-KC20

User Guide

Contents

1 Introduction **2**
1.1 Alcohol and Its Effects on the Human Body 2
1.2 Alcohol Is a Drug 2
1.3 How Does Alcohol Affect My Body? 2
1.4 How Does Alcohol Affect My Driving? 3

2 Dose-Specific Effects of Alcohol **4**

3 Sensor Cleaning **4**

4 Preparation **5**
4.1 Installing the Battery 5
4.2 Components Diagram 5

5 Operation **5**
5.1 Error Notifications 6

6 Specifications **7**

7 Precautions **7**

8 Verification and Calibration **7**
8.1 Calibration Instructions 8

9 Limits of Breath Alcohol Testing **8**

10 Service and Support **8**
10.1 Manufacturer 9

1. Introduction

This device is intended to measure alcohol in human breath. Measurements obtained by this device are used in the diagnosis of alcohol intoxication.

BACtrack Go is an alcohol screening device, used for the detection of alcohol in the breath. The BACtrack provides a digital result, indicating the approximate BAC (Blood Alcohol Content) of the test subject. The BACtrack is powered by one AAA battery and is very easy to use. Simply turn on the unit and begin testing in a matter of seconds. However, before you begin testing, there are a few things you should know. Please take a moment to read your user guide in its entirety and familiarize yourself with the BACtrack.

1.1. Alcohol and Its Effects on the Human Body

Alcohol is absorbed from the mouth, throat, stomach and intestines into the bloodstream.

Alcohol ingested by the human body can be detected in the breath because of its presence in the bloodstream.

Alcohol cannot be digested, and it cannot be chemically changed in the bloodstream. As the blood flows through the lungs, alcohol in the bloodstream moves across the membranes of the lung's air sacs (alveoli) into the air. The concentration of the alcohol in the alveolar air is directly related to the concentration of the alcohol in the blood. As the alveolar air is exhaled, the alcohol in it can be detected by the breath alcohol testing device.

The alcohol concentration in the breath is related to that in the blood, and because of this, an individual's BAC can be determined by measuring alcohol in the breath. The ratio of breath alcohol to blood alcohol is generally estimated to be 2,100:1. Therefore, 2,100 milliliters (ml) of alveolar air will contain approximately the same amount of alcohol as 1 ml of blood.

1.2. Alcohol Is a Drug

Alcohol is the chemical (ethanol or ethyl alcohol) resulting from the fermentation of grapes or grain. Alcohol is absorbed directly into your bloodstream. From your bloodstream, the fluids in your body tissues absorb the alcohol. Your brain is made up of a large concentration of fluids and will absorb a large amount of any alcohol you drink. Your liver eliminates the absorbed alcohol at its own rate and nothing you eat or drink can speed up the process. Your liver does this by oxidizing the alcohol (converting it into water and carbon dioxide). Coffee, food or any other "cure" will not sober you up. Only time will do that.

1.3. How Does Alcohol Affect My Body?

Alcohol is a depressant. It has a relaxing effect on the muscles in your body. The muscles of your eyes relax and lose focus. Your eyesight will become fuzzy and you

may experience double vision. Brain activity is slowed. Your judgment, reflexes, and coordination are all negatively affected.

Some vision impairments that occur when you have been drinking include:

- narrowing of your field of vision;
- reduction in your depth perception;
- decreased ability to see in darkness;
- increased sensitivity to glare and a longer time for your eyes to readjust from the glare.

Some mental (brain) impairment can occur when you have been drinking, including:

- reduced awareness of danger;
- becoming overly confident and reckless;
- difficulty in making decisions;
- reduction in balance;
- slowed reflexes;
- impaired judgment.

1.4. How Does Alcohol Affect My Driving?

We can't emphasize enough how dangerous it is for you to drink and drive. Your vision and your brain are the most important factors in driving safely. If you jeopardize either by drinking, and then you drive, you are likely to be involved in a serious or fatal collision.

Source: California Department of Motor Vehicles, USA (`http://www.dmv.ca.gov`).

2. Dose-Specific Effects of Alcohol*

BAC	Dose-specific effects[1]
0.02–0.03 %	Slight euphoria. Loss of shyness. Depressant effects are not apparent. Impairment possible in some individuals.
0.04–0.06 %	Feeling of well-being, relaxation, lower inhibitions and sensation of warmth. Euphoria. Some minor impairment of reasoning and memory. **Lowering of caution. Driving skills may be impaired at this level of intoxication.**
0.07–0.09 %	Slight impairment of balance, speech, vision, reaction time, and hearing. Euphoria. Judgment and self-control are reduced. Caution, reason and memory are impaired. **Driving skills are always impaired at this level of intoxication and higher.**
0.10–0.12 %	Significant impairment of motor coordination and loss of good judgment. Speech may be slurred. Balance, vision, reaction time and hearing will be impaired. Euphoria.
0.13–0.15 %	Gross motor impairment and lack of physical control. Blurred vision and major loss of balance. Euphoria is reduced and dysphoria (anxiety, restlessness) is beginning to appear.
0.16–0.20 %	Dysphoria predominates. Nausea may appear.
0.25 %	Needs assistance in walking. Total mental confusion. Dysphoria with nausea and some vomiting.
0.30 %	Loss of consciousness.
0.40 % and up	Onset of coma. Possible death due to respiratory arrest.**

[1] Bailey, William J., *Drug Use in American Society*, 3rd ed., Minneapolis: Burgess, 1993.

* The effects of alcohol intoxication are greatly influenced by individual variations among users. Some users will be intoxicated at a much lower BAC than shown above.

** Death can occur at a lower BAC in some individuals.

The generally accepted legal standard for alcohol intoxication in the United States is 0.08 %. However, your driving skills can be impaired at a level of 0.04 % or lower. It is never safe to drink any amount of alcohol and drive.

3. Sensor Cleaning

When the breathalyzer is used for the first time, or if the unit has not been used for an extended period of time (3 weeks or more), airborne pollutants may form on the sensor surface. To maximize accuracy, clean the sensor by performing 2 tests, allowing it to run a complete cycle each time, and disregard the results. Whether or not alcohol is present on the breath for the sensor cleaning doesn't matter. After


running through the complete testing process twice, cleansing is complete, and you can test again normally.

4. Preparation

4.1. Installing the Battery

Install one AAA alkaline battery in the battery compartment.

4.2. Components Diagram

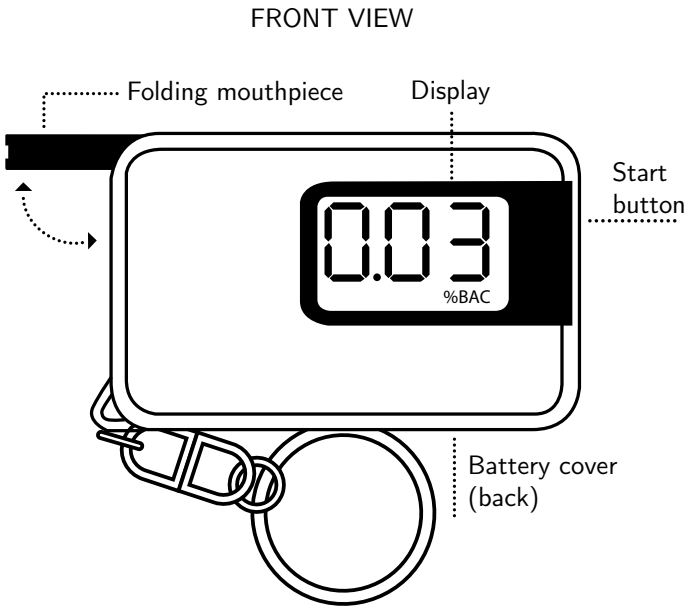5. Operation

IMPORTANT!

Test subjects should wait 20 minutes after eating, drinking, or smoking before blowing into the BACtrack breathalyzer. **Failure to wait 20 minutes can provide inaccurate test results and damage the sensor.**


Step 1: (Optional) Attach a mouthpiece cover to the folding mouthpiece.

Step 2: Press the start button. The BACtrack will begin a countdown.

Step 3: When the countdown approaches 2, begin to inhale a deep breath.

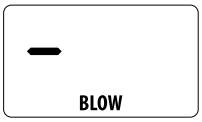Step 4: When the countdown reaches zero, the BACtrack will beep. Blow through the mouthpiece continuously for five seconds until the beep sound stops.

Step 5: Once the sensor has analyzed the breath sample, the estimated BAC value is displayed.


Note: For most accurate results, wait one minute between tests. When testing high alcohol levels, wait 3–5 minutes between tests.

5.1. Error Notifications


Display shows ***FL***
The FL notification can appear in three cases:

- The user did not blow immediately after the countdown.
- The user did not blow continuously for five seconds.
- The user may still have alcohol in their saliva, and the measured alcohol concentration is uneven. Users should wait 20 minutes after their last drink to make sure alcohol in the mouth is fully absorbed and a proper breath sample can be taken.

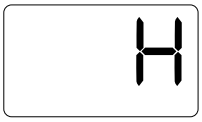Display shows ***H***
Reading exceeds upper limit of 0.40 % BAC. Make sure the user has waited 20 minutes after consuming alcohol before blowing.

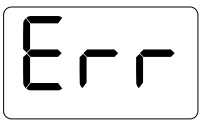Display shows ***Err***
The sensor is damaged. Contact support for assistance (see section 10).

Low battery icon
When the low battery icon appears, replace the battery with a new AAA alkaline battery.

6. Specifications

Dimensions	1.42 × 0.56 × 2.39 in (36.16 × 14.14 × 60.77 mm)
Weight	2 oz (56 g)
Sensor technology	MicroCheck™
Detection range	0.00–0.40 % BAC
Power supply	One AAA alkaline battery (included)
Battery life	50–75 tests
Warm-up time	11–13 seconds
Response time	5 seconds
Operating temperature	50–104 °F (10–40 °C)

7. Precautions

1. Wait at least 20 minutes after drinking, eating or smoking before testing. Not observing this waiting period can cause inaccurate readings and damage the sensor.
2. Do not blow smoke, food, or liquids into the BACtrack, as this will damage the sensor.
3. Do not test in areas with strong winds, smoke, or in an area where a large amount of alcohol is being consumed.
4. Send your tester in for periodic calibration service as required (see section 8).
5. The BACtrack breathalyzer is designed to be used in a temperature range of 50–104 °F (10–40 °C).
6. Avoid testing in the presence of any substances that contain methyl alcohol, isopropyl alcohol or acetone. These substances may interfere with the results of the test.
7. You cannot use the results of this product in court.
8. Do not use the BACtrack as a tool to determine whether you should operate a motor vehicle or equipment, or perform any other dangerous act.
9. Do not drink and drive. Always have a designated driver when alcohol is being consumed.

8. Verification and Calibration

The BACtrack is calibrated during manufacture using advanced alcohol simulation equipment. Known alcohol concentrations are passed through the sensor to set

baseline values for testing. During testing, the unit compares users' breath samples to these baseline values.

The accuracy of breath alcohol testers can fluctuate after six to twelve months of normal use, depending on the number of tests performed and operating conditions. If this product is providing inconsistent test results, is not providing any test results, or provides unusually high or low test results, the product should be recalibrated immediately by a trained professional.

The product should be recalibrated at least every six to twelve months. If the product is used every day, it may need to be recalibrated as frequently as every month.

8.1. Calibration Instructions

To have your BACtrack recalibrated, contact your local distributor or the store where you bought it. The manufacturer's information on recalibration is available at `http://support.bactrack.com` or by phone on +1 415-693-9756 (USA).

During recalibration, a professional technician will reset the sensor baseline values and also perform a full diagnostic check to ensure the accuracy and reliability of the product.

9. Limits of Breath Alcohol Testing

This breath alcohol tester is intended to measure alcohol in human breath. Readings obtained by this device are used in the diagnosis of alcohol intoxication.

This product cannot determine whether, or to what extent, a test subject's mental or physical functioning or judgment may be impaired by blood alcohol, including whether the test subject is intoxicated under any definition of that word. Any act or failure to act based on a reading from this product shall be at the test subject's or user's own risk. KHN Solutions assumes no responsibility for consequences to or of test subjects who use this product and are later shown to have been under the influence of alcohol or to have had their judgment or any mental or bodily function impaired by alcohol.

Correlation between breath alcohol content and blood alcohol concentration depends on many variables, including environmental factors (such as air quality, wind, humidity, temperature, etc.) and health conditions of the test subject. A low BAC reading does not mean that the test subject's physical or mental performance or judgment (including reaction times) can respond to an emergency. The concentration of alcohol in the blood of a test subject cannot be exactly determined by using a breath alcohol screening product.

10. Service and Support

If you have any questions, please contact your local store or distributor, or visit `https://www.bactrack.com`.

If the breathalyzer needs to be repaired or recalibrated, contact your local distributor for more information.

10.1. Manufacturer

KHN Solutions Inc
300 Broadway, Suite 26
San Francisco, CA 94133, USA
`https://www.bactrack.com`

Distributor
Sunnysoft s.r.o.
Kovanecka 2390/1a
190 00, Praha 9 - Liben
www.sunnysoft.cz

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